

# Exam Room Equipment

Product range for Examination and Procedure Rooms includes:

Heine, Halux and MaggyLamp Exam Lights, Magnifier Light, X-Ray Viewers, MediRoo Lock Box, Seca Scales, Seca Height Measures, Exam Tables Fixed and Electric, Gynae Exam Tables, Surgeon Stools / Chairs, Step Stools and Privacy Screens



Suppliers of quality products to the Healthcare Industry est 1889.  
Australian Owned & Operated.

Web [dbhealth.com.au](http://dbhealth.com.au)

Email [customerservice@denyer.com.au](mailto:customerservice@denyer.com.au)

Freecall: 1800 451 446

ABN: 37 000 459 174

## Heine EL10 LED Examination Light

The Heine EL 10 Examination Light distinguishes itself through an absolutely homogenous field of light. The light intensity can be regulated continuously between 6,500 and 45,000 Lux. The spot size can be adjusted in a diameter between 8.5 and 14cm to suit any situation at a working distance of 30cm.

- Controls with ergonomic grip allow for intuitive and fast adjustment of brightness and spot size. The compact illumination head (Ø approx. 60mm) allows almost coaxial illumination.
- Bright and homogeneous. Perfect illumination in every examination situation.
- 10 Watt high performance LED with authentic colour rendering
- Easy and efficient cleaning and disinfection. Easy and intuitive adjustments via coiled goose neck.
- Electrical control system monitors all functions and guarantees safety.

Description	Product Code
EXAM LIGHT - HEINE EL10 WHITE LED W/ WALL MT (EACH)	HJ0827001
EXAM LIGHT - HEINE EL10 WHITE LED W/ CLAMP MT (EACH)	HJ0827002
EXAM LIGHT - HEINE EL10 WHITE LED W/ MOBILE STAND (EACH)	HJ0827003



HJ0827003



Adjusting spot size on EL10

## Heine EL3 LED Examination Light

- The EL3 LED provides bright shadow free illumination precisely at the point of interest. At a working distance of 30cm > 42,000 lux / >350 lumens.
- The LED illumination provides optimal light intensity, homogeneity and colour rendering for the most accurate diagnosis.
- Compact Illumination Head (diameter 6cm) allows for convenient, precise positioning, particularly in difficult application situations.
- Cool Touch closed head design makes disinfection safe and easy.
- The lightweight head stable goose neck allow easy and intuitive adjustment to keep the light in the required position.

Description	Product Code
EXAM LIGHT - HEINE EL3 WHITE LED W/ WALL MT (EACH)	HJ0827011
EXAM LIGHT - HEINE EL3 WHITE LED W/ CLAMP MT (EACH)	HJ0827013
EXAM LIGHT - HEINE EL3 WHITE LED W/ MOBILE STAND (EACH)	HJ0827014



HJ0827011

## Halux LED Exam Light

The Halux LED exam lights are made in Switzerland by Derungs. These quality lamps have been a proven performer in surgeries, clinics and hospitals around Australia. The Halux LED is an ergonomically designed, easy to use light which gives 40,000 lux at 0.5 metre.

- 40,000 Lux at 0.5 metres, 180mm Ø beam width and 7 Powerful LEDs.
- 4400 K (cool white) ideal for recognition of veins or wound treatment.
- Easy to clean smooth edges.
- 950mm reach with Spring-balanced arm system for precise light positioning.
- On/off rocker Switch on arm.



Description	Product Code
EXAM LIGHT - HALUX LED 20-1 W/ MOBILE STAND (EACH)	RBML508M
EXAM LIGHT - HALUX LED 20-1 W/ WALL MOUNT (EACH)	RBML508W

## Maggylamp Examination Light—Halogen and LED

Maggylamp Halogen Examination Lights are an entry level range of lamps that offers a cost effective solution for surgery examination lighting.

- 700mm gooseneck.
- Sure grip handle with on/off switch.
- Warm white colour 2300K colour.
- 2 years warranty.
- Wall or mobile mounted (View all options below).
- 4,000 hour bulb life.
- 45,000 lux at 30cm.
- 35 watt, 24° beam.



RBML308

Maggylamp LED Examination Light with its flexible gooseneck arm and specifically selected colour corrected LED's offers flexibility, correct light output and colour required for surgery examinations.

- 40,000 hour lamp life.
- Cool white light, 5500K colour.
- Low energy consumption, <5 watt.
- 2 years warranty.
- Low heat output.
- 3 x 1 watt LED giving 16,000 lux at 30cm.
- Concise beam angle of 30°.
- Wall or mobile mounted.
- 700mm gooseneck.



RBML308LED

Description	Product Code
EXAM LIGHT - MAGGYLAMP HALOGEN (EACH)	RBML308
EXAM LIGHT - MAGGYLAMP LED (EACH)	RBML308LED
OPTION—WALL BRACKET RIGHT ANGLE FOR MAGGYVUE & MAGGYLAMP(EACH)	RBML307
OPTION—MOBILE STAND FOR MAGGYVUE & MAGGYLAMP #95CM (EACH)	RBML130

## Maggyvue Opal Magnifier Lamp

The Maggyvue Magnifier Lamp is a quality German designed lamp. It is very popular with the medical, electronics and inspection industries where sharp definition is required.

- Quality German optics.
- Circular lens, Ø 120mm, 2 x magnification.
- Optional clip on lens to increase magnification #ML151.
- Spring balanced joints, easy to adjust.
- Supplied with bench clamp and diffuser.
- 22 watt daylight colour fluorescent tube.



RBML150

Description	Product Code
MAGNIFIER LAMP—MAGGYVUE OPAL (EACH)	RBML150
REPLACEMENT TUBE—22W FOR ML150 (EACH)	RBML235
OPTION—BENCH CLAMP FOR MAGGYVUE #0-65MM (EACH)	RBML110
OPTION—WALL BRACKET RIGHT ANGLE FOR MAGGYVUE & MAGGYLAMP (EACH)	RBML307
OPTION—MOBILE STAND FOR MAGGYVUE AND MAGGYLAMP #95CM (EACH)	RBML130
OPTION—CLIP ON LENS FOR MAGGYVUE #3X MAGNIFICATION (EACH)	RBML151

## Dalcross Slimline LCD X-Ray Viewers

- Surface Mounted Single, Double and Triple Bay Slimline LCD X-Ray Viewers
- Stylish design with 2 key-hole slots in rear for wall mounting
- Supplied with plug lead that exits at bottom right of light box
- High brightness and whiteness with even, flicker-free illumination
- Aluminium housing with good heat dissipation
- Positive grip holds X-ray film secure
- Maintenance free. No globes to change



DA1601SL



DA1600SL

Description	Product Code
X-RAY VIEWER - SINGLE SLIM LINE LCD #47 X 54.5X 2.4CM (EACH)	DA1600SL
X-RAY VIEWER - DOUBLE SLIM LINE LCD #85 X54.5 X2.4CM (EACH)	DA1601SL
X-RAY VIEWER - TRIPLE SLIM LINE LCD #123 X54.5 X2.4CM (EACH)	DA1602SL

## Dalcross Halogen X-Ray Viewer

- Surface Mounted Single, Double and Triple Bay Halogen X-ray Viewers
- 2 Key-hole slots in rear for wall mounting
- Single and double can be upright on base, for use on desk or bench. Triple Bay viewer is Wall Mounted only
- Supplied with plug lead that exits at bottom right of light box
- Excellent light diffusion and uniform illumination is achieved from: circular fluorescent tube, continuous opal perspex viewing panel and internal white reflector surface
- Stainless steel strip contains grip rollers to hold X-rays and Each individually switched bay is internally divided



DA1601

Description	Product Code
X-RAY VIEWER - SINGLE #37.7X50X10.5CM (EACH)	DA1600
X-RAY VIEWER - DOUBLE #75X50X10.5CM (EACH)	DA1601
X-RAY VIEWER - TRIPLE #112.5X50X10.5CM (EACH)	DA1602



DA1602

Also available in Flush Mount and on a mobile stand.

Available in 4 in a line, 6 in a line, 8 in a line, spinal X-Ray, in 2 tiers models—2 x 2, 3 x 3 and 4 x 4.

## MediRoo Wall Lock Box

MEDIroo Lock Box are designed for securing medications, narcotics, supplies, personal belongings, and other valuables at the point-of-care or anywhere secure storage is needed.

- Constructed of durable 20-gauge powder-coated steel - can be cleaned with bleach and common surface disinfectants.
- Can be easily mounted to a wall.
- Doors open at a 180 degree angle and do not protrude into the room and door locks automatically when door is closed.

Small—352 x 113 x 206mm

Large—597 x 113 x 206mm



Description	Product Code
MEDIROO LOCK BOX—SMALL KEYED LOCK WHITE (EACH)	CR53005-81
MEDIROO LOCK BOX—LARGE KEYED LOCK WHITE (EACH)	CR53006-81



## Seca Baby Electronic Scales

**Seca 334 Baby Scale** is designed as a reliable portable baby scale with a integrating weighing tray larger and is more comfortable than other models.

**Seca 354 Baby Scale** is designed as a dual function portable baby / child scale. At a touch of a button, the tray can be detached quickly and easily from the base leaving you with robust flat scales for children weighing up to 20 kg. As the two elements are securely connected with one another, it is absolutely safe to lift and transport the scales.

**Seca 374 Baby Scale** is designed as a dual function portable baby / child scale with extra large weighing tray. Seca 374 is also compatible with 360 wireless system - allowing measurements to be sent wireless to the wireless printer or to analytics 115 on your PC.



PN334

Model	Seca 334	Seca 354	Seca 374
Capacity & Graduation	20kg - 5gm up to 10kg and 10gm for 10kg plus	20kg - 10gm up to 10kg and 20gm for 10kg plus	20kg - 5gm up to 10kg and 10gm for 10kg plus
Functions	Automatic switch-off, TARE, Auto-HOLD, BMI Function	Breast Milk Intake, TARE, Auto-HOLD and BMI Function	TARE, Breast Milk Intake Function, auto-HOLD, auto-CLEAR, damping, automatic switch-off, SEND/PRINT
Weight & Dimensions	2.8kg - 63.8 x 30 x 10.5cm	2.3kg - 55.2 x 33.2 x 15.6cm	3.6kg - 62 x 35.8 x 19cm
Power	Batteries	Batteries	Batteries
Options	400 Power Adaptor, 232 Measuring Rod and 428 Carrying Case	413 Carrying Case	400 Power Adaptor, 233 Measuring Rod, 418 Head and Foot Positioners and 360 compatible

Description	Product Code
SCALE - DIGITAL BABY SCALE 20KG (EACH)	PN334
SCALE - DIGITAL BABY SCALE 20KG (EACH)	PN354
SCALE - DIGITAL BABY SCALE WITH X-LARGE TRAY 20KG (EACH)	PN374
SWITCHMODE POWER SUPPLY UNIT—SECA ADAPTOR FOR 703, 634, 869, 334 AND 374 (EACH)	PN400
OPTION - SECA MEASURING ROD FOR BABY SCALES (EACH)	PN232
CARRY BAG FOR BABY SCALES #334 (EACH)	PN428
TRANSPORT CASE—SECA FOR BABY SCALE 354 (EACH)	PN413
MEASURING ROD—SECA FOR BABY SCALE 374 #35-80CM (EACH)	PN233
HEAD AND FOOT POSITIONER—SECA FOR BABY SCALE 374 (EACH)	PN418



PN354



PN374

## Seca Baby Measuring Mat

- Seca 210 Measuring Mat is kind to the skin, washable and foldable.
- The light weight, space-saving solution for easy, precise measurement of the length of babies and toddlers while lying down.
- The fixed head piece and the easy to slide foot positioner make it simple to use.
- Measuring range: 10 to 99cm
- Weight: 580g

Description	Product Code
MEASURING MAT - SECA BABY 10 - 99CM (EACH)	PN210



PN210

## Seca Platform and Wheelchair Electronic Scales

The range of platform scales below have the same scale but with different features such as handrails and ramps for wheelchairs.

- Seca 634 Electronic Platform Bariatric Scale - Cable remote display with 2.5m cable can be held, wall mounted or placed on desktop. Operates from batteries but is available with optional 447 Power Adaptor.
- Seca 676 Electronic Wheelchair Scale - includes 2 integrated ranges and handrail with connected display. Operates from power adaptor.



PN634



PN676

Model	634	284
Capacity & Graduation	360kg - 100gm graduation	360kg - 100gm graduation
Functions	TARE, pre-TARE., auto-HOLD, BMI, kgs/lbs switch-over, damp-ing, Automatic Switch-off	TARE, Pre-TARE, auto-HOLD, BMI, auto-CLEAR, kgs/lbs switch-over, dampening , auto switch-off, SEND/PRINT, Transport Castors, 2 x Integrated Ramps.
Weight & Dimensions	14.4kg - 64 x 60 x 5.3cm	38kg—92 x 115 x 112cm
Power Supply	Batteries	Power Adapter
Options	472 Stand for cable remote display, and power adaptor 447.	N/A

Description	Product Code
PLATFORM SCALE - SECA ELECTRONIC 360KG (EACH)	PN634
STAND FOR CABLE REMOTE DISPLAY SECA (EACH)	PN472
SWITCHMODE POWER SUPPLY UNIT—SECA ADAPTOR FOR 703, 634, 869, 334 AND 374 (EACH)	PN400
PLATFORM SCALE - SECA ELECTRONIC 360KG W/ RAIL & RAMPS	PN676

## Seca Chair Scales

The Chair Scales are designed to be able maintain the patient's comfort while providing a reliable chair scale.

- Seca 954 has a high capacity of 300kg and is powered by power adaptor and rechargeable battery.
- The 954 is compatible with the 360 wireless system.
- Additionally the 954 has some extra functions including BMI and SEND/ PRINT.
- Seca 952 has a capacity of 200kg and is powered by batteries.
- Both chairs include swivelling armrests, secure wheel brakes and foldaway footrests.



PN954



PN952

Description	Product Code
SCALE - SECA DIGITAL CHAIR C/W BMI CALC. 300KG (EACH)	PN954
SCALE - SECA DIGITAL CHAIR 200KG (EACH)	PN952

## Seca Electrical Column Scales

- Seca 703 Column Scale is extremely robust, precise and stable with a capacity of 300 kg,
- a large but low platform and the
- Integrated automatic BMI function,
- these column scales have the prerequisites for the diagnosis and therapy of heavyweight patients.

<b>Model</b>	<b>703</b>
Capacity & Graduation	300kg - 50gm up to 150kg than 100gm graduation
Functions	TARE, pre-TARE, mother/child function, HOLD, auto-HOLD, BMI, auto-CLEAR, damping, Automatic Switch-off, SEND/PRINT
Weight & Dimensions	15.2kg - 36 x 52 x 93cm
Power Supply	Batteries
Options	447 Power Adaptor, and 220 Height Measuring Rod



Description	Product Code
SCALE - SECA DIGITAL COLUMN 300KG #WIRELESS (EACH)	PN703
HEIGHT MEASURE - SECA TELESCOPIC ROD FOR SCALES (EACH)	PN220
SWITCHMODE POWER SUPPLY UNIT—SECA ADAPTOR (EACH)	PN400

## Seca Wireless Measuring Station

Wireless measuring station Seca 284 measures and weighs in 1 step.

- Scale with capacity of 300kg (50g graduation) and a height measuring range of 30-220cm (1mm graduation).
- Multi-function touch display enables to the features on the measuring station, a second display is on the height measure.
- Both Scale and Height Measure are RS232 compatible and can wireless send result via 360 system.

<b>Model</b>	<b>284</b>
Capacity & Graduation	300kg—50gm graduation
Functions	TARE, Pre-TARE, Mother/Child Function, HOLD, auto-HOLD, auto-BMI, CLEAR, auto-CLEAR, auto weighing range switch-over, dampening and SEND/PRINT. Height Measure Functions: HOLD, auto switch-off, user-defined zero setting and SEND/PRINT.
Weight & Dimensions	16kg—43.3 x 46.6 x 239.4cm
Power Supply	Power Adapter
Options	N/A



Description	Product Code
MEASURING STATION - SECA WIRELESS SCALE & STADIOMETER #300KG / 220CM (EACH)	PN284
MOBILE MEASURING SYSTEM—SECA (EACH)	PN254

## Seca Mobile Measuring Station

Seca 254 Mobile Measuring Systems integrated:

- Electronic Flat Scales with 200kg capacity (graduation 50gm up to 150kg and than 10gm),
- stadiometer with height range of 20-205cm,
- Adaptor Element to connect scale and stadiometer
- Carry case.
- Designed for mobility this set enables weighing and height measurement for medical professionals.

Description	Product Code
MOBILE MEASURING SYSTEM—SECA (EACH)	PN254



## Seca Floor Scales

**Seca 813 Floor Electronic Scale** has designed with a wide, low platform and non-slip dimpled mat. Making the Seca 813 easy and safe to mount also for older or handicapped patients. Due to its high capacity of 200kg and stable integrated steel frame it is also ideal for weighing very heavy patients.

**Seca 762 Floor Mechanical Scale** 150kg has designed with a large, low-level, non-slip platform and large well defined dial designed for everyday use in medical clinics.

**Sea 803 Floor Electronic Scale** has designed a general use scale (150kg) with a quality rubber surface and large LCD display.

**Seca 876 Floor Electronic Scale** is a robust, light weight scale with mother/child function and is regularly used for mobile use.



PN813



PN762

Model	762	803	876	813
Capacity & Graduation	150kg - 500gm	150kg - 100gm	250kg - 100gm to 150kg than 200gm	200kg - 100gm
Functions	Reset to Zero	auto-HOLD, tap-on automatic switch-on and switch-off	Mother-Child Function, tap-on automatic switch-on and switch-off	Auto-HOLD, tap-on automatic switch-on and switch-off
Power	Mechanical	Batteries	Batteries	Batteries

Description	Product Code
SCALE - SECA MECHANICAL PRECISION DIAL 150KG (EACH)	PN762
SCALE - SECA DIGITAL FLOOR 150KG #BLACK (EACH)	PN803B
SCALE - SECA DIGITAL FLOOR 150KG #WHITE (EACH)	PN803W
SCALE - SECA DIGITAL FLOOR WIDE / LOW 200KG (EACH)	PN813
SCALE - SECA DIGITAL FLOOR 250KG #MOBILE USE (EACH)	PN876



PN876



PN803

## Seca Height Measures and Stadiometers

**Seca 206 Wall Mounted Height Tape Measure** is the space saving solution with roll-up mechanism. This tape measure is designed for uncomplicated installation with just screw required to used.

**Seca 222 Wall-mounted Height Measure** has been specially developed for wall mounting, this measuring rod is aligned to work exactly to the millimetre from head to toe.




PN206

**Seca 213 Stadiometer** is designed for mobile use and generally used to measure children and teenagers without the restriction of having a height measurement wall mounted.

**Seca 274 Wireless Stadiometer** is designed for mobile use. This system is designed to be compatible with Seca 360 wireless making recording measurement simple.

Description	Product Code
HEIGHT MEASURE - SECA 2.2M TAPE WALL MT (EACH)	PN206
HEIGHT MEASURE - SECA TELESCOPIC ROD W/HEEL #WALL MT (EACH)	PN222
STADIOMETER - SECA PORTABLE HEIGHT MEASURE 20-205CM (EACH)	PN213
STADIOMETER—SECA ELECTRONIC PORTABLE HEIGHT MEASURE	PN274

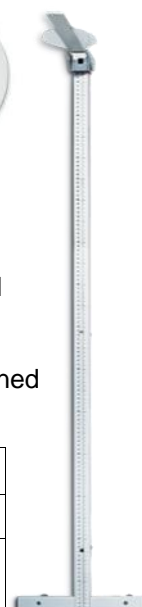

PN222



PN213



PN274



## Seca 360 Integration with EMR

A large number of scales and height measures are now capable of transmitted measurements that can be feed into your EMR.

### SECA 456—USB WIRELESS ADAPTOR

- With the seca 456 you equip your PC to receive wireless data from seca scales and measuring systems.
- You can eliminate annoying cable connections in the examination rooms and you can store measurements in the seca database, feed them into an EMR system, and provide other parameters in a PDF file (EMR-integrated). The labor-intensive process of documenting examination findings is turned into a speedy and paperless matter.



PN456

### SECA 101—SECA EMR FLASH 101

- Free software to connect Seca 360 wireless products to EMR systems.

### SECA 466—WIRELESS PRINTER

Seca 466 is a wireless printer for the seca 360° range of products to receive, analyse and print measurements on thermal paper or labels.

- Network capable via seca 360 wireless technology
- Analysis and graphic formatting of measurements
- Comparison of height and weight
- Children's results expressed in percentile



Scales that can be integrated in Seca 360 include:  
374, 284, 703, 274, 676, 954.

Description	Product Code
SECA 360 USB ADAPTER (EACH)	PN456
SECA 360 WIRELESS PRINTER (EACH)	PN466

## Fixed Height Examination Couch

The Fixed Height Exam Couch are cost effective exam couches essential for consulting and treatment rooms. Combining comfort, stability, strength and aesthetic design ensure that these exam tables remain an essential part of your surgery.

Features of the DA1000 2 Section Exam Couch are:

- Sturdy rectangular steel frame and easy to assemble—flat packed (4 bolts)
- Adjustable backrest with Vinyl covered top available in a Grey or Black.
- Available in Grey Powder Coat
- Dimensions: 185 x 61 x 80cm (High) and Max User Weight 175kg



DA1000

Features of the DA1004 2 Section Heavy Duty Exam Couch are:

- Sturdy rectangular steel frame and easy to assemble—flat packed (boltless)
- Adjustable backrest with Vinyl covered top available in a Grey or Black.
- Available in Grey Powder Coat
- Dimensions: 183 x 61 x 78cm (High) and Max User Weight 250kg



DA1004

Description	Product Code
EXAM COUCH - HEAVY DUTY P/C #250KG SWL (EACH)	DA1004
EXAM COUCH - STANDARD P/C #175KG SWL (EACH)	DA1000

## High-Low Examination Couch

The High-Low Examination Couches provide a cost effective alternative for medical professionals on a budget.

- Safe working load - 200kg
- Sturdy and reliable heavy duty construction.
- Quiet electronic motor controlled by foot-switch operation for effortless positioning.
- Retractable castors available on all models.
- Height adjustment from 46 to 91cm.
- Dimensions (approx.): 195.5 x 71cm.



Retractable Castors



DA1242E



DA1243M



DA1243E

Description	Product Code
EXAM COUCH - ELECTRIC HI-LO P/C #2 SECTION W/ MANUAL HEADREST (EACH)	DA1242M
EXAM COUCH - ELECTRIC HI-LO P/C #2 SECTION W/ ELEC HEADREST (EACH)	DA1242E
EXAM COUCH - ELECTRIC HI-LO P/C #3 SECTION (EACH)	DA1243M
EXAM COUCH - ELECTRIC HI-LO P/C #3 SECTION #3 MOTORS (EACH)	DA1243E

## Multipurpose Hi-Lo Exam Chair

All electric couch: one motor for HiLo, one motor for seat tilt and trendelenburg. The seat tilt provides extra leg elevation for podiatry. The leg rest goes to near vertical for comfortable sitting.

- Footswitch controlled electric motors for Hi-Lo lift and trendelenburg tilt functions.
- Three section upholstered top with gas strut assisted backrest and legrest
- Overall dimensions: 630mmW x 1820mmL
- Height: 500mm to 950mm
- Head tilt angles: 85° to -15° and Leg tilt angles: 15° and -85°
- 225kg Safe working load
- Strong, rigid frame with dual lifting points.
- Quality upholstery finish in blue vinyl
- Flip back armrests
- Maintenance free, no need for regular lubrication or adjustment

Bariatric Version also available—350kg SWL—DA1251



Description	Product Code
EXAM COUCH - ELECTRIC HI-LO P/C MULTIPURPOSE CHAIR (EACH)	DA1248

## Transforma Gynae Exam Couch

Transforms from a Exam HiLo couch to specialist Gynae Table without removing or adding anything. Suitable for use as a dual purpose outpatients couch.

- The mid section lifts up to become a backrest in gynae position
- 3 Motors for HiLo, tilt and midsection backrest
- Power height adjustment 500mm to 950mm
- Electronic tilt for colposcopy or trendelenberg to 15°
- Head & mid section tilt +85°
- Gas-strut end section headrest
- Adjustable padded footpads
- Stainless steel debris tray
- Retractable castors
- Black or blue vinyl upholstery
- Overall Dimensions: 1830mmL x 700mmW
- Headrest: 640mmL x 630mmW
- Backrest (mid section): 790mmL x 700mmW
- Seat section: 400mmL x 700mmW
- Safe working load: 200kg



Description	Product Code
EXAM COUCH - ELECTRIC HI-LO P/C TRANSFORMA GYNAE (EACH)	DA1261

## Nina Gynae Exam Couch

The ultimate table for OB/GYN examination and procedures. No other product can match the functionality and quality finish of the Nina.

- Positions: Hi-Lo, pelvic tilt, trendelenberg, forward tilt
- Power height adjustment 480mm to 880mm
- Three section contoured seamless upholstery
- Removable leg extension stows-away behind the backrest
- Adjustable lithotomy footpads
- Dual column lift and motorized backrest
- Electronic footswitch and hand controls
- Debris tray
- Adjustable head cushion
- Locking wheels
- Overall Dimensions: 1840mmL x 700mmW
- Backrest: 900mmL x 700mmW
- Seat section: 480mmL x 700mmW
- Leg extension: 430mmL
- Safe working load: 200kg



Description	Product Code
EXAM COUCH - ELECTRIC HI-LO P/C TNINA GYNAE (EACH)	DA1259G

## Gas Lift Stools with and without Fixed Back

The Gas Lift Stools have been designed for medical professional use.

- 37cm diameter seat with vinyl upholstery
- Height adjustment from 53 to 72cm via gas lift lever located under seat
- Option: Choice of either chair castors or chair glide.
- Choice of vinyl upholstery colour - standard colour White, Blue or Black



OX10-3232



OX10-3234

Description	Product Code
STOOL - OXFORD BACKLESS GAS LIFT #NYLON BASE (EACH)	OX10-3230
STOOL - OXFORD BACKLESS GAS LIFT #ALLOY BASE (EACH)	OX10-3234
STOOL - OXFORD FIXED BACK W/ GASLIFT #NYLON BASE (EACH)	OX10-3231
STOOL - OXFORD FIXED BACK W/ GASLIFT #ALLOY BASE (EACH)	OX10-3236
STOOL - OXFORD W/ ADJUSTABLE FOOT RING #BACKLESS #NYLON BASE (EACH)	OX10-3232
STOOL - OXFORD W/ ADJUSTABLE FOOT RING #BACKLESS #ALLOY BASE (EACH)	OX10-3235
STOOL - OXFORD GAS LIFT W/ FIXED BACK REST & FOOTRING #NYLON BASE	OX10-3233
STOOL - OXFORD GAS LIFT W/ FIXED BACK REST & FOOTRING #ALLOY BASE	OX10-3237



OX10-3235

## Deluxe Gas Lift Chair

The Deluxe Gas Lift Chair has been designed for general hospital use.

- Comfortable ergonomically shaped seat and backrest.
- Height adjustment (53 to 72cm).
- Backrest adjustment (height/angle - 10 to 16°).
- Adjustable foot ring and cast alloy base has 65cm diameter.
- Choice of either chair castors or chair glide.
- Choice of vinyl upholstery colour

Description	Product Code
CHAIR - OXFORD DELUXE GAS LIFT W/ FOOTRING #NYLON BASE (EACH)	OX10-3213
CHAIR - OXFORD DELUXE GAS LIFT W/ FOOTRING #ALLOY BASE (EACH)	OX10-3214



OX10-3213

## Apex Surgeon Stool with Hand & Foot Control

- Height adjustment 550mm to 750mm and Dual controls for seat height adjustment by hand or foot
- Triangular shaped seat allows better freedom of movement of users legs, increasing maneuverability
- Height, angle adjustable backrest and Seat tilt adjustment - Seat Diameter 420mm
- "Brake unloaded" castors prevent the stool from moving when unoccupied but become freewheeling when weight is on the chair - Five star alloy base
- Premium vinyl with easy clean plastic covers
- Vinyl upholstery colours: blue
- Safe working load: 120kg

Description	Product Code
STOOL - APEXSURGEON W/ HAND & FOOT CONTROL #ALLOY BASE (EACH)	DA1408





## Bambach Saddle Chair with Foot Control

Designed to reproduce the upright standing position, whilst sitting down. It's specially moulded shape holds your pelvis in this upright neutral position, which allows your spine to fall into perfect alignment.

- Height adjustment 606mm to 730mm
- Foot control for hands-free height adjustment
- Strong and stable five star moulded base with 50mm rubber castors
- Height and tilt adjustable backrest with contoured lumbar support
- Safe working load: 140kg



Description	Product Code
STOOL - DALCROSS BAMBACH SADDLE SEAT W/ BACKREST (EACH)	DA1403
STOOL - DALCROSS BAMBACH SADDLE SEAT (EACH)	DA1402

## Foot Operated Surgeon's Stool and Blood Collection Chair

- Foot operated adjustable seat height (550 to 740mm)
- Comfortable 350mm diameter stool pad
- Castors standard: 50mm twin



The Blood Collection Chair has been designed specifically for pathology centres and maybe used in number of applications.

- This chair features:
- Foot step and Wide angled armrests (20cm Max).
- Choice of vinyl upholstery colour - standard colour White, Blue or Black
- Seat Width: 44cm and Seat Height: 62cm



Description	Product Code
SURGEON'S STOOL—OXFORD FOOT OPERATED P/C (EACH)	OX10-1570
CHAIR - OXFORD BLOOD COLLECTION W/ FOOTREST (EACH)	OX09-3018

## Step Stool

The Single and Double Step Stools are designed to meet demands for Exam Room such as getting on or off exam couches.

- These Step Stools are made from tubular powder coated frame with non-slip aluminium tread surface.
- The Surgeons step stools are made from a stainless steel frame with a non-slip aluminium top.
- These step stools can be easily customised to meet the requirements of individual users.



Description	Product Code
STEP STOOL - OXFORD SINGLE S/S #40X25X19CM (EACH)	OX10-2166
STEP STOOL - OXFORD DOUBLE P/C #58X42X40CM (EACH)	OX10-2160
FOOT STOOL - OXFORD SURGEON'S SINGLE #40X30X20CM (EACH)	OX10-3850



## Step Stools with Handrails

- The Single and Double Step Stools with Handrail are designed to meet demands for patients with mobility issues.
- These Step Stools are made from tubular powder coated frame with non-slip aluminium tread surface.

Description	Product Code
STEP STOOL - OXFORD SINGLE P/C W/ HANDRAIL #57X30X24-100CM (EACH)	OX14-2145
STEP STOOL - OXFORD DOUBLE P/C W/ HANDRAIL #70X60X40-135CM(EACH)	OX14-2164
STEP STOOL - OXFORD DOUBLE P/C W/ TWIN HANDRAIL #70X60X40-135CM(EACH)	OX14-2163



OX14-2145

OX14-2164

OX14-2163

## Privacy Screens

- Frames with optional curtains
- Easy removable washable curtains
- Available in a range of powder coated finishes with curtain colours to compliment
- Standard colours: light blue curtains and white birch powder coating
- Castors standard: 50mm twin
- 2 Fold Privacy Screen—1800 x 480 x 1640mm high
- 3 Fold Privacy Screen—2700 x 480 x 1640mm high



Description	Product Code
PRIVACY SCREEN - OXFORD 3 FOLD W/ CURTAIN P/C (EACH)	OX13-2201
PRIVACY SCREEN - OXFORD 2 FOLD W/ CURTAIN P/C (EACH)	OX13-2217